ProCuity Disposable Seizure / Fall Risk Side Rail Pads

INSTALLATION GUIDE | ITEM 116-7237



1 Open package

Open the package. Remove the pads and covers. Allow foam to expand 2-5 minutes.



Packaged materials




Unpacked foam pads

ALLOW 2-5 MINUTES

2 Match each pad

Ensure rails are in up position. Match each cover with foam to its matching rail. The cover's stamped side should always face caregiver, and left/right refers to patient's anatomical sides.

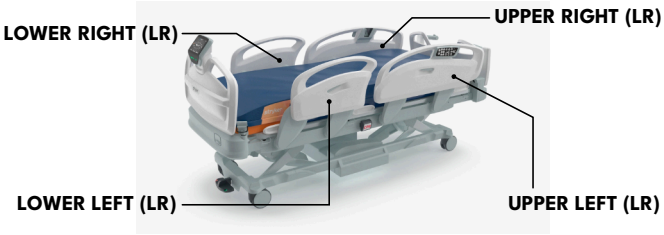


3 Open the straps

Gently separate the fastening straps:
3 fastening straps per upper pad;
2 fastening straps per lower pad.



4 Slip the fitted cover with foam over rail

Slip the cover and foam fully over its matching rail, with stamped side facing caregiver and foam side towards patient.

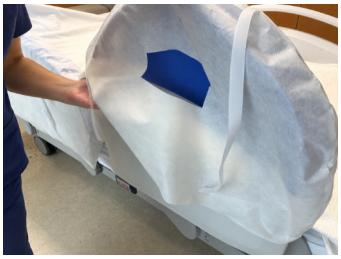

5 Guide the bottom straps

Bring top/front straps toward you. Hold the top of the pad. Gently feed bottom/back straps down between the mattress and rail.



6 Fasten and align

Match and secure each corresponding pair of straps. Align openings with rail panels.



Upper pad: 3 fastening straps. Lower pad: 2 fastening straps.

FINAL CHECK

- ☐ Pads in place
- ☐ All straps secure
- ☐ Rail panels accessible

Global Medical Foam, Inc. **globalmedfoam.com** **(419) 884-9354**